

Dynamic Measurements of Greenhouse Gas Respirations Caused by Changing Oxygen Levels

Derek Fleck^{1*}, Nabil Saad¹

1. Picarro, Inc., Santa Clara, CA, United States. (*Contact: dfleck@picarro.com)

Abstract

The necessity for constant monitoring of greenhouse gases (GHGs) is clearly evident now more than ever. Moreover, interpreting and understanding the processes that dictate the production and consumption of these gases will allow for proper management of GHGs in order to mitigate its detrimental climate effects. Presence of oxygen, or lack of it, is the driving force for determining pathways within biochemical redox reactions. Experiments to find correlations between oxygen and greenhouse gases have helped us understand photosynthesis, denitrification and beyond. Within the past few years measurements of O₂ and nitrous oxide have been used over a wide ranging array of disciplines; from studying avenues for redox chemistry to characterizing gas profiles in sputum of cystic fibrosis patients. We present a full analysis solution, based on cavity ring-down spectroscopy, for simultaneous measurements of N₂O, CO₂, CH₄, H₂O, NH₃, and O₂ concentrations in soil flux, in order to better understand dynamics of ecological and biogeochemical processes.

The stability and high temporal resolution of the five-species CRDS analyzer, coupled with a continuous high-precision O₂ measurement (1-σ <200ppm) produces a complete picture of biogeochemical processes, for which a multitude of additional research experiments can be conceived. Adding another dimension to explore to help determine the rate at which these greenhouse gases are produced or consumed, allows scientists to further address fundamental scientific questions. Data is presented showing precision, drift and limitations of the O₂ sensor measurement as well as the validity of spectroscopic corrections with the CRDS analyzer caused by changing O₂. Experimental data is also presented to explore correlations of gas flux rates of O₂, N₂O, CO₂ and CH₄ in a plant/soil ecosystem due to differing soil and atmospheric conditions at varying timescales from minutes to days.

Instrumentation

Picarro G2508 CRDS Gas Concentration Analyzer



- ~ Continuous gas flow at ~240 sccm
- ~ Continuous concentration measurement of N₂O, CH₄, CO₂, NH₃ and H₂O in air with ppb precision
- ~ Data rate: Measures all 5 species in < 6 seconds

Completely Integrated Electro-optical O₂ Sensor - PreSens

- ~ Pressure and temperature stabilized

Peripheral Sensors for Environmental Properties

- ~ Soil O₂ concentration - Apogee S-210
- ~ Ground Temperature - Apogee
- ~ Ambient Temperature - Campbell T109
- ~ Solar Radiation - Licor Li190R

Methodology


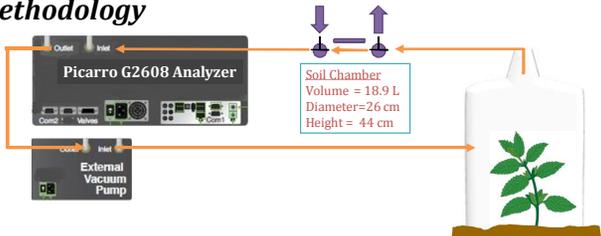


Figure 1: Back of the Picarro analyzer and vacuum pump showing gas connections to the chamber two three way valves (purple) selected to vent chamber automatically

- ~ A closed-loop system using 1/4" Swagelok Teflon tubing and Swagelok fittings recirculates air between the chamber and the Picarro analyzer
- ~ The chamber was closed for 60 min for each measurement
- ~ Chamber vented for minimum of 1 hr (later data 2 hrs) to bring chamber closer to ambient concentrations (though full equilibrium never reached)

Analysis


Analyzer Stability and Drift – Allan Standard Deviation

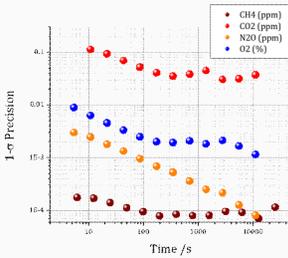


- ~ Long Term Stability of instrument due to pressure and temperature stabilization allows low long term drift
- ~ More importantly short term drift is completely negligible allowing for no calibration to be done throughout experiment
- ~ Maximum Precision for O₂ is around 5 minutes averaging and can reach 20ppm

O₂ Correction and Validity


- ~ Changing O₂ concentration changes shape of absorption feature used for concentration
- ~ Changes of 5% can cause an influence of more than 0.5% of the concentration measured
- ~ Changes are small but influence can be seen when data is averaged
- ~ Current corrections are valid over a range of O₂ levels, from nearly 5% to 30% total O₂

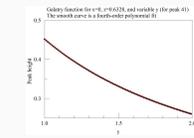

**** Background:** These deviations stem from spectral features line shape changing and apply to all absorption spectrometers. Figure to the Left shows simulation of changing peak height because of y-parameter changing. This can be caused by things such as Changing O₂:N₂ ratios.

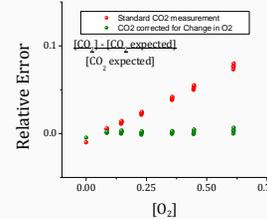


Raw Concentration Measurements for Flux Determination

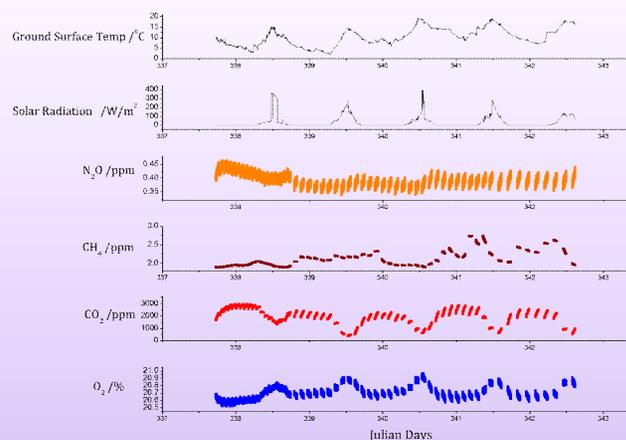


Raw Concentration measurements show diurnal cycles that the CO₂ and O₂ concentrations follow. Also show that chambers were not fully vented to ambient concentrations, but since the fluxes of CO₂ and O₂ are large, can still be easily seen

Results

Diurnal Cycles of Plant Ecosystem: Production and Consumption

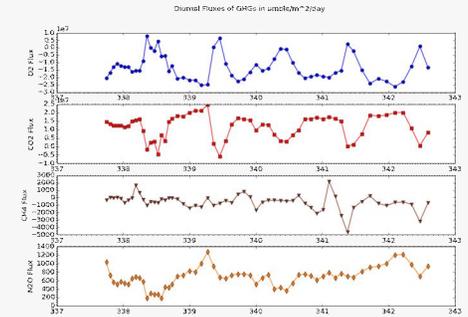

Above Figure shows the fluxes for all gases measured over the period of 5 days. It can clearly be seen that there is a strong correlation to day and night time, as well as between CO₂ and O₂. All fluxes are calculated in μmole/day/m²

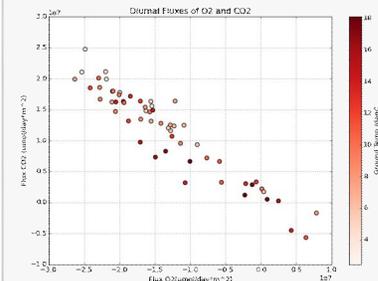

The correlation between CO₂ fluxes and O₂ fluxes is shown to the left. This small single plant ecosystem only has net positive fluxes of O₂ at peak solar radiation, but lower ground temperatures, so microbial activity is not at its height consuming O₂. This can explain the large ground temperatures in the middle where there is a balance of consumption from soil with production of from the plant.

This produced a fitted slope of - 0.77 CO₂ flux to O₂ flux.

Conclusions and Recommendations

In this study, we have demonstrated a integrated sensor to a proven instruments in the CRDS Picarro G2508, This new implementation allows automatic corrections to be made within the analyzer due to widely varying O₂:N₂ ratios which can be caused by microbial consumption or plant production. Additionally This O₂ sensor has allowed us to see the changing background matrix to the normal GHGs monitored in soil respiration experiments. We have found the microbial consumption of O₂ can be so large that it can even offset the O₂ production due to a mint plant. Likewise CO₂ production due to soil microbes significantly outweighs the consumption of CO₂ from the mint plants, except at the sunniest part of the day in winter at 450 W/m². It should also be noted that, no useful data can be seen in diurnal patterns of the soil oxygen content. soil O₂ measurements using the Apogee, which could be due to long response and poor digitization of signal.